# LinguSil

Gebruiksaanweisung

Instructions for use

Mode d'emploi

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**D**  
DENTAURUM

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## LinguSil

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# LinguSil

## Instructions for use

Thank you for having chosen this quality product from Dentaureum.

In order to use this product at its best for you and your patients, it is important to study and follow these directions for use carefully.

The written instructions cannot cover all eventualities during operation. For this reason our Hotline is available to answer any other questions and ideas that may arise.

Our products are regularly modernized to the latest technical development. For this reason we advise to always check the most up-to-date operating instructions supplied for any possible changes, even if the product is already being used frequently.

### Manufacturer

Dentaureum GmbH & Co. KG

Turnstraße 31 · 75228 Ispringen · Germany

### LinguSil, REF 164-000-00

### Polyvinyl silicate for direct application

#### Areas of application

- Transfer matrix for brackets (indirect bonding)
- Fixing of x-ray balls during implant planning
- Front-, occlusal stamps


#### Material characteristics

- Application time: 45 seconds
- Working time -  
including oral  
working time: 3 minutes
- Shore A Hardness: 72
- Product color: transparent

## Instructions for use

This product is designed for single use only!

After having inserted the cartridge into the **LinguSil Injector** (REF 164-001-00), remove the protective cap by turning.

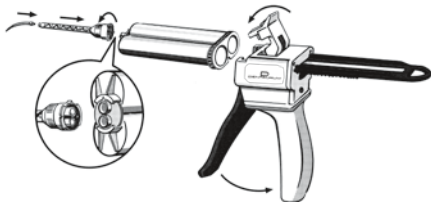LinguSil Injector, REF 164-001-00

LinguSil, REF 164-000-00

**Before use always ensure that the cartridge opening is not blocked.** Press handle so that a pea sized amount emerges and until both components come out evenly. Assemble the mixing tip and fix in place by turning to the right. The dosage and mixing will occur automatically when the cartridge is pressed. Replace the protective cap after use (hygiene).

## Caution

Set silicon material is permanent, for this reason, avoid contact with clothes.

The setting reaction of all polyvinyl silicate materials are affected by retroactive solutions with iron chloride or aluminum chloride bases as well as methacrylates, latex gloves and sulfur compounds. It is therefore essential to avoid contact with these materials.

## Temperature

The specific working and oral working times are based upon room temperatures (23°C/73.4°F) and the temperature within the mouth (32°C/89.6°F). Higher temperatures reduce and lower temperatures increase these working times.

**Storage: 18 - 25°C / 64.4 - 77°F**

**Shelf life: 24 months; when stored in a cool environment approx. 3-6 months longer**

## Disposal

This product can be disposed of in the domestic refuse.

**The Dentaureum orthodontic hotline is available for you should you have any further questions concerning our products.**

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